



Chair Scale

USER MANUAL **MS5410**




Please keep the instruction manual at hand and follow instruction for use.

CONTENTS

Explanation of Graphic Symbols on Label/Packaging ...	3
I. Safety Notes	6
A. General Information.....	6
B. EMC Guidance and Manufacturer's Declaration.....	10
II. Installation.....	14
A. Adjusting footrests	14
B. Adjusting wheel height	16
C. Inserting Batteries	17
D. Using Adapter	20
III. Indicator	21
A. Indicator and Key Functions	21
B. Display layout.....	22
IV. Using Device	23
A. Basic Operation	23
B. Hold	23
C. BMI.....	24
D. Tare	24
E. Pre-Tare	25
F. Print	29
V. Device Setup	30
A. Setting Time & Date	30
B. Device Setup	31
VI. Setup USB Connection to PC	33
VII. Wireless Connection.....	36
VIII. Troubleshooting.....	37
Error Messages.....	38
IX. Product Specifications	39
A. Device Information.....	39
B. Power Adapter Standards.....	40
X. Declaration of Conformity	43

Explanation of Graphic Symbols on Label/Packaging

Text/Symbol	Meaning
	Caution. Failure to follow instructions may result in adverse events
	Separate collection for waste of electrical and electronic equipment, in accordance with Directive 2002/96/EC. Do not dispose of device with everyday waste
	Name and address of device manufacturer
	Carefully read user manual before installation and usage, and follow instructions for use
	Medical electrical device, Type B applied part
	Medical electrical device, Type BF applied part
	Device catalogue number
	Name and address of authorized representative in the European Union (applicable for medical devices only)
	Device is a medical device. Text indicates device category type
	Manufacturer's batch or lot number for device
	Device's serial number
	Device's model number
	Device's Unique Device Identifier (01) GTIN (11) Production date (21) Serial number
	Value in mass units. This is the difference between two consecutive display values, used to classify and verify a scale
	Device conforms to (EU) 2017/745 Regulation on Medical Devices. Four digit number is identifier for medical device Notified Body

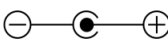	Device claims conformity with essential regulations set by the relevant European Union legislation
	Device complies with EU directives M: Conformity label in compliance with Directive 2014/31/EU for non-automatic weighing instruments 17: Year in which conformity verification was performed and the CE label was applied. (ex: 17=2017) 0122: Identifier for metrology Notified Body
	Device is a Class III scale in compliance with Directive 2014/31/EU
	Name and address of entity importing device
	Name and address of entity responsible for translating Information For Use
	Event counter recording how many times protected-mode parameter adjustments relating to metrology-related characteristics (ex: gravity value) have been conducted for device
	Device has received Taiwan NCC approval
	Device conforms to U.S. Federal Communications Commission regulations
	Device complies with all UK applicable product legislation
	Device's polarity of power.
	The Chair Scale is for weight measurement only. Please do not push the Chair Scale while the patient is sitting on it.
	The footrests are provided for the subject to rest their feet upon, for accurate weight measurement. Do not stand on the footrests!

"In case of differences, icon on device itself takes precedence"

Copyright Notice

Charder Electronic Co., Ltd.

No.103, Guozhong Rd., Dali Dist., Taichung City 41262 Taiwan

Tel: +886-4-2406 3766

Fax: +886-4-2406 5612

Website: www.chardermedical.com

E-mail: info_cec@charder.com.tw

Copyright© Charder Electronic Co., Ltd. All rights reserved.

This user manual is protected by international copyright law.

All content is licensed, and usage is subject to written authorization from Charder Electronic Co., Ltd. (hereinafter Charder) Charder is not liable for any damage caused by a failure to adhere to requirements stated in this manual. Charder reserves the right to correct misprints in the manual without prior notice, and modify the exterior of the device for quality purposes without customer consent.



Charder Electronic Co., Ltd.
No. 103, Guozhong Rd., Dali Dist.,
Taichung City, 41262 Taiwan

I. Safety Notes

A. General Information

Thank you for choosing this Charder Medical device. It is designed to be easy and straightforward to operate, but if you encounter any problems not addressed in this manual, please contact your local Charder service partner.

Before beginning operation of the device, please read this user manual carefully, and keep it in a safe place for reference. It contains important instructions regarding installation, proper usage, and maintenance.

Intended Purpose

This medical device is designed to be used in accordance with national regulations, to measure weight within specifications, for weight-related usage by professionals.

Clinical Benefit

Measurement results can be used by professionals to diagnose (and monitor) weight-related issues.

Intended medical indications/contraindications

Measurement: patient's body weight. No known contraindications to measurement of body weight.

Intended patient profile

- (a) Age: no restrictions
- (b) Weight: no restrictions within device weight capacity
- (c) Patient Conditions: require measurement of body weight.
Unable to stand independently without support.

Intended user profile

- (a) At least 20 years old
- (b) Minimum knowledge:
 - To be able to read at a high-school level and understand Arabic numerals (e.g. 1, 2, 3, 4...)
 - Basic hygiene knowledge
 - Trained in device's operation
 - Read the instruction manual

(c) Language

- Able to read the language of instruction manual and on-screen instructions

(d) Qualifications

- No special certifications or qualifications required

Residual Risk Evaluation

- (a) All foreseeable risks have been evaluated and considered acceptable. Generally speaking, the most likely risk caused by incorrect usage of the device is less accurate measurement (or inability to use device to acquire measurement), which does not pose imminent physical risk to patient or user.
- (b) Benefit-risk ratio is considered acceptable. Chair scales are an important option for measuring patients. Usage of device is unlikely to result in harm to user or patient.

General Handling

- Ensure all parts are properly locked and tightened before operating the device.
- Measurement accuracy requires the subject's feet, back, and head to be straightly aligned. Please note that height can vary throughout the day
- **CAUTION:** Do not use next to equipment that may cause electromagnetic or other types of interference.

Safety Instructions

Before putting device into use, please read this user manual carefully. It contains important instructions for installation, usage, and maintenance of device.

The manufacturer shall not be liable for damages caused by failure to heed the following instructions:

- The device has an expected service life of 5 years when correctly handled, serviced, and periodically inspected in accordance with manufacturer's instructions.
- Improper installation will render the warranty null and void.
- Observe permissible ambient temperatures for use

Cleaning

- Device surface should be cleaned using alcohol-based wipes.

Maintenance

- Please contact your local Charder distributor for regular maintenance and calibration, regular checking of accuracy is recommended; frequency to be determined by level of use and state of device.

Warranty/Liability

- The period of warranty shall be eighteen(18) months, beginning on the date of purchase. Please retain your receipt as proof of purchase.
- No responsibility shall be accepted for damage caused through any of the following reasons: unsuitable or improper storage or use, incorrect installation or commissioning by the owner or third parties, natural wear and tear, changes or modifications, incorrect or negligent handling, chemical, electrochemical, or electrical interference, unless damage is attributable to negligence on the part of Charder.
- This device does not contain any user-maintained parts. All maintenance, technical inspections, and repairs should be conducted by an authorized Charder service partner, using original Charder accessories and spare parts. Charder is not liable for any damages arising from improper maintenance or usage. Dismantlement of the device will void the warranty.



Warning

Measurements for physically disabled persons.

- Physically disabled persons should not attempt to take measurements alone, but instead should have their caretakers assist them in using the device.
- Footrest can only be used when subject is sitting in chair. To avoid injury, subject should refrain from standing on footrest, as device may tip over if used incorrectly.



Incident Reporting

- Any serious incident that has occurred in relation to the device should be reported to the manufacturer, EU representative (if device is used in EU member state), and competent authority of user/subject's member state.

B. EMC Guidance and Manufacturer's Declaration

Guidance and manufacturer's declaration-electromagnetic emissions		
The product is intended for use in the electromagnetic environment specified below. The customer or the user of the product should assure that it is used in such an environment.		
Emission test	Compliance	Electromagnetic environment-guidance
RF emissions CISPR 11	Group 1	The product is suitable for use in all establishments other than domestic and those directly connected to a low voltage power supply network which supplies buildings used for domestic purposes.
RF emissions CISPR 11	Class A	
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations /flicker emissions IEC 61000-3-3	Compliance	

Guidance and manufacturer's declaration-electromagnetic immunity

The product is intended for use in the electromagnetic environment specified below. The customer or the user of the product should assure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment-guidance
Electrostatic discharge(ESD) IEC 61000-4-2	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%
Electrical fast transient/burst IEC 61000-4-4	±2kV for power supply lines	± 2kV for power supply lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	±1kV line(s) to line(s) ±2kV line(s) to earth	±1kV line(s) to line(s) ± 2kV line(s) to earth	Mains power quality should be that of a typical commercial or hospital environment.
Voltage Dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0% UT for 0,5 cycle 0% UT for 1 cycle 70% UT(30% dip in UT) for 25cycles 0% UT for 5 s	0% UT for 0,5 cycle 0% UT for 1 cycle 70% UT(30% dip in UT) for 25cycles 0% UT for 5 s	Mains power quality should be that of a typical commercial or hospital environment. If the user of the product requires continued operation during power mains interruptions, it is recommended that the product be powered from an uninterruptible power supply or a battery.
Power frequency(50, 60 Hz) magnetic field IEC 61000-4-8	30 A/m	30 A/m	The product power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.

NOTE UT is the a.c. mains voltage prior to application of the test level.

Guidance and manufacturer's declaration-electromagnetic immunity

The product is intended for use in the electromagnetic environment specified below.

The customer or the user of the product should assure that is used in such and environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment-guidance
Conducted RF IEC 61000-4-6	3 Vrms 150 KHz to 80 MHz <u>6 V in ISM bands between 0.15 MHz and 80 MHz</u> <u>80 % AM at 1 kHz</u>	3 Vrms 150 KHz to 80 MHz <u>6 V in ISM bands between 0.15 MHz and 80 MHz</u> <u>80 % AM at 1 kHz</u>	Portable and mobile RF communications equipment should be used no closer to any part of the product including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.
Radiated RF IEC 61000-4-3	3 V/m <u>80MHz to 2,7 GHz</u>	3 V/m <u>80MHz to 2,7 GHz</u>	Recommended separation distance: $d = 1,2 \sqrt{P}$ $d = 1,2 \sqrt{P}$ 80MHz to 800 MHz $d = 2,3 \sqrt{P}$ 800MHz to 2,7GHz Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres (m).
			Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, ^a should be less than the compliance level in each frequency range. ^b
			Interference may occur in the vicinity of equipment marked with the following symbol:
			

NOTE1 At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the product is used exceeds the applicable RF compliance level above, the product should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the product.

b Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3V/m.

**Recommended separation distance between
portable and mobile RF communications equipment and the product**

The product is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the product can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the product as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter m		
	150 kHz to 80 MHz $d = 1,2\sqrt{P}$	80 MHz to 800 MHz $d = 1,2\sqrt{P}$	800 MHz to 2,7 GHz $d = 2,3\sqrt{P}$
0,01	0,12	0,12	0,23
0,1	0,38	0,38	0,73
1	1,2	1,2	2,3
10	3,8	3,8	7,3
100	12	12	23

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in metres (m) can be estimated using the equation applicable to the frequency of the transmitter, where p is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.


NOTE1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.


II. Installation

A. Adjusting footrests


1. Remove bolt to adjust length of footrest

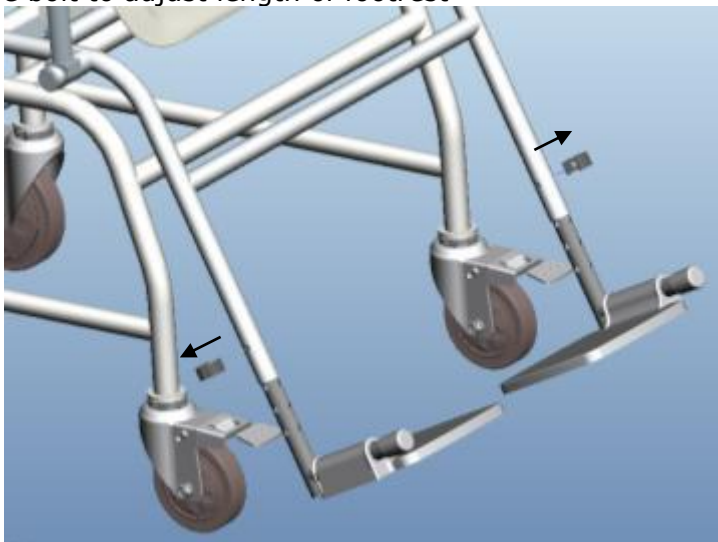

2. Adjust footrest height according patient leg length

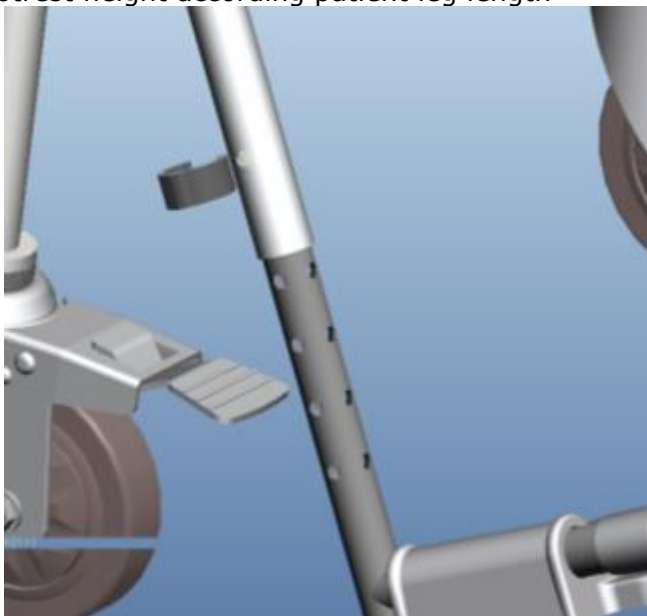

3. After completing adjustment, insert bolts to secure footrest.

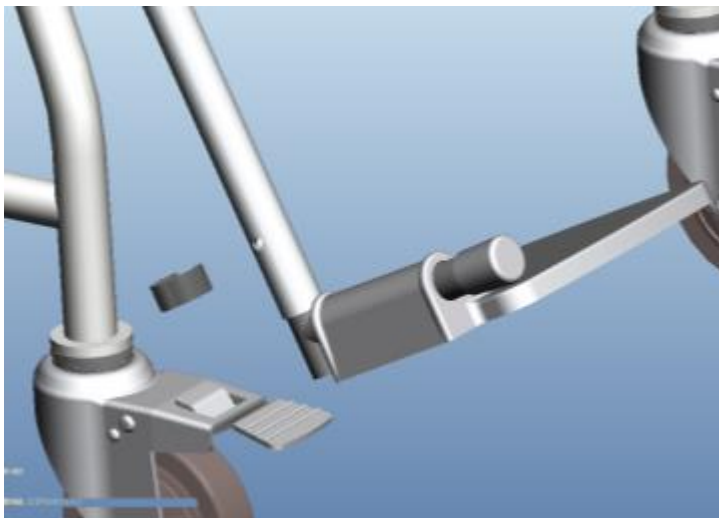

B. Adjusting wheel height

1. Place device on flat surface, apply wheel brakes



Brakes


2. To tighten wheel castor, loosen counter nut slightly. After loosening, turn wheel frame clockwise to tighten.



Counter nut

Wheel frame


3. Adjust wheel height until air bubble on level indicator is level

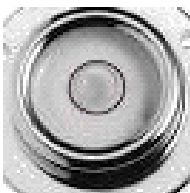



NOTE: Be careful not to lose wheels during adjustment

C. Inserting Batteries

1. Open battery housing cover

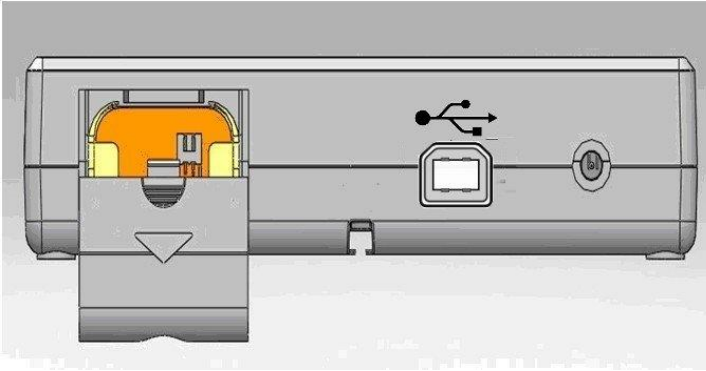


2. Accessing batteries

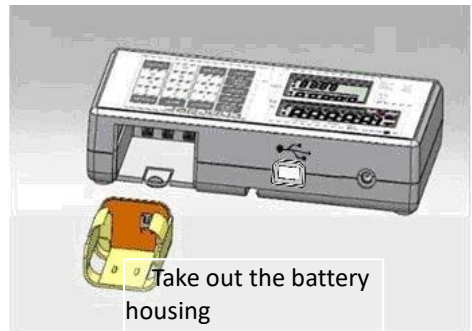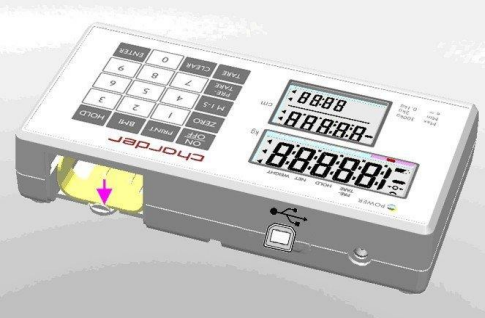


Take out the battery housing


3. Use either rechargeable battery pack, or AA batteries

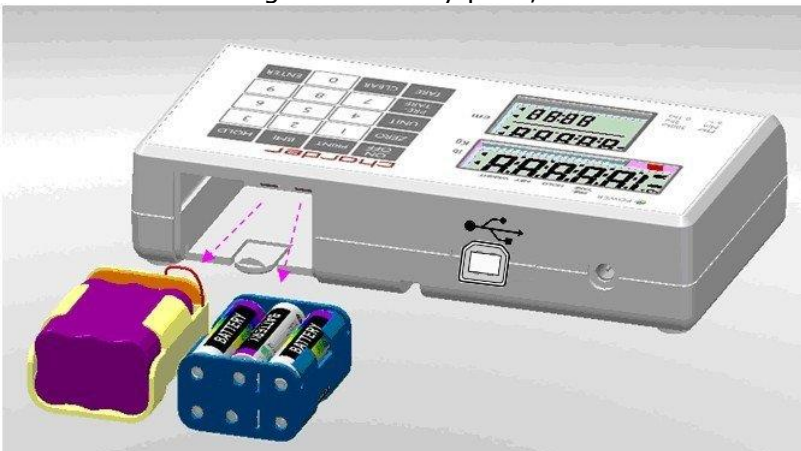


4. Ensure batteries are installed into the housing correctly

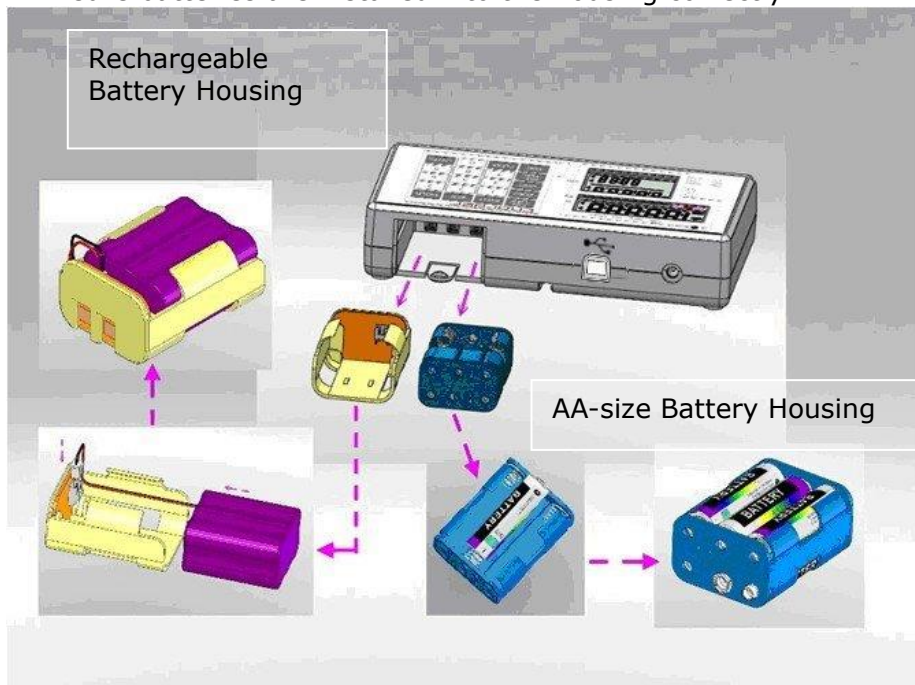

5. Install the battery housing into the compartment, and make sure the right side of housing pin is facing towards inside of the connecting position

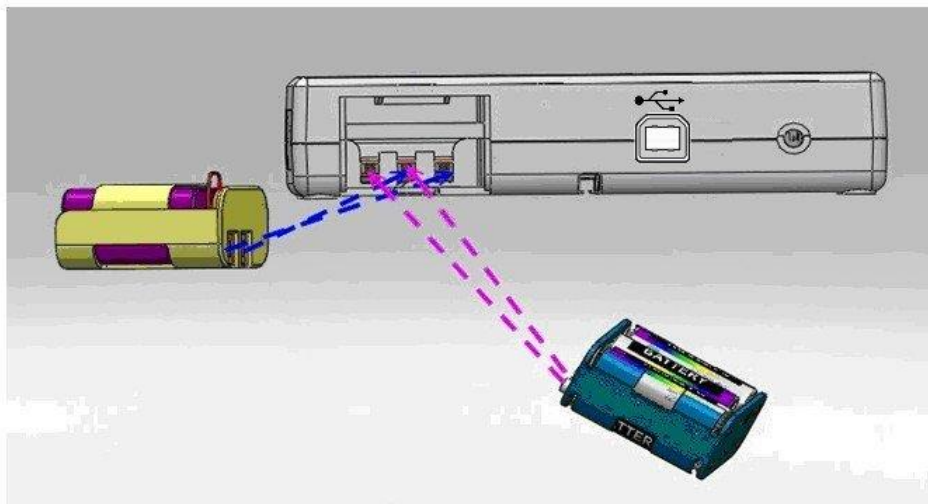

6. Slide back the cover to close the battery housing compartment. Turn on power to confirm that battery is correctly installed.

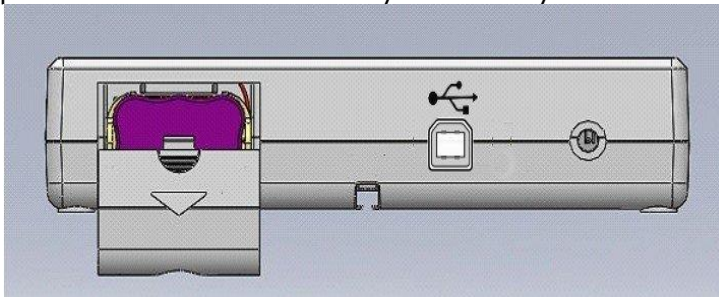


Using Rechargeable Battery (optional)

The rechargeable battery should be recharged at least once every 3 months, regardless of if the device has been used. Battery can be charged by plugging device's exclusive adapter into AC Connector Port.

After a long period in storage (e.g. >3 months), the battery should run a full cycle (charge/discharge) to allow it to restore full capacity.


Ensure rechargeable battery housing is installed and inserted properly into the compartment.

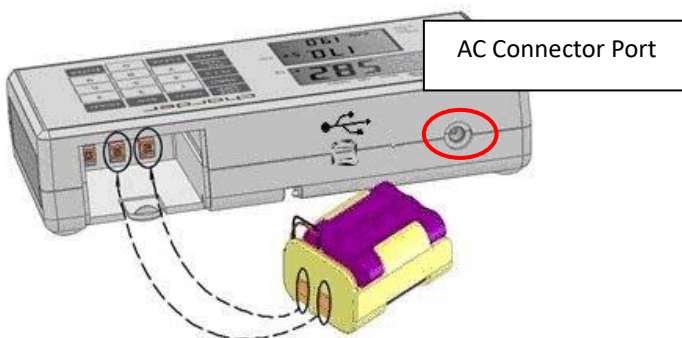

If . Lo prompt displays on the LCD, please charge battery promptly to avoid battery damage.

D. Using Adapter

1. Connect adapter to indicator before connecting to mains power supply
2. Disconnect adapter from mains power supply before unplugging adapter pin from indicator.

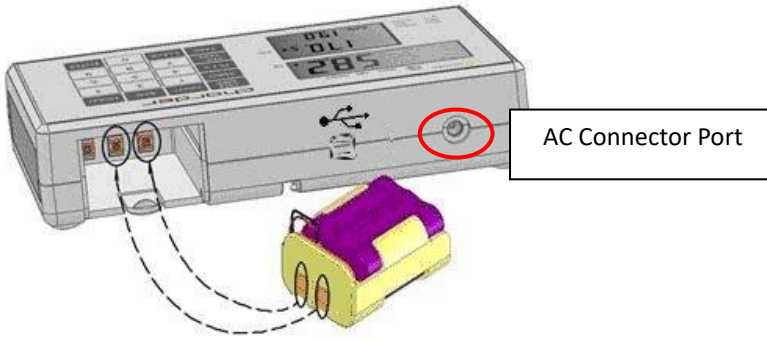

III. Indicator

A. Indicator and Key Functions




(Wireless functionality optional)

Key Function

ON/OFF	Power on or power off.
ZERO	Reset display to 0.0 kg display. Press and hold for 3 seconds to enter device settings.
M1-5	Save pre-tare values (up to 5)
PRE-TARE	Pre-tare the known weight of an object (ex: chair) before beginning measurement.
TARE	Allows user to deduct weight from reading after measurement
PRINT	When printer or PC is connected to the scale, press this key to print results
BMI	Calculation of Body Mass Index
HOLD	Determine stable weighing value - used when weight is unstable. Press and hold for 3 seconds to enter time setting.
0-9	For entering digits.
CLEAR	Clear incorrect data input.
ENTER	Confirm input.

B. Display layout

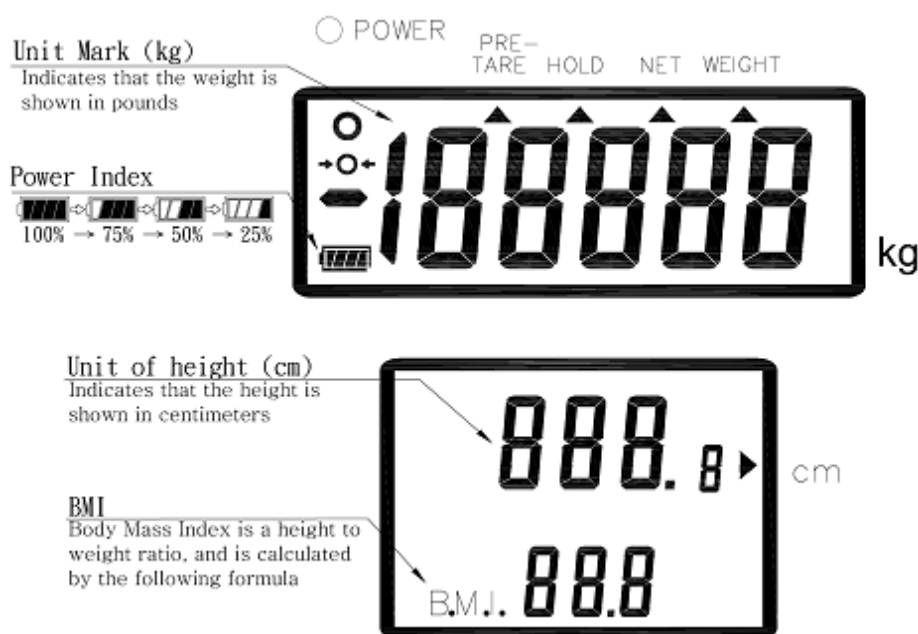

Definitions

Stable symbol: Indicate that weight is stable.

Zero symbol: Weight is at zero

Negative weight: Weight under zero.

Low battery: Battery needs to be charged or replaced.

IV. Using Device

A. Basic Operation

Switch on the device using  key. The device will automatically perform self-calibration, displaying software version.

Once "0.00 kg" appears on indicator, device is ready for measurement.

Note: If "0.00 kg" does not display on indicator, press  key to zero the device.

Guide subject to sit on chair. Make sure subject's feet are off the ground, and properly placed upon footrests. After the weight has stabilized, the "stable" symbol will appear on indicator.

Note: If subject's weight exceeds scale capacity, indicator will display "Err" prompt due to overload.

B. Hold

The hold function determines average weight, designed to be used if subject's weight will not stabilize (ex: an active child).

Note: If fluctuation is too severe, average weight determination will be difficult and hold may not function correctly.

1. Switch on the device normally.

2. Press the  key. "HOLD" will be displayed on the indicator.

3. Guide subject to sit on chair.

4. After a few seconds, the average weight will be displayed on the indicator. This weight will be locked - at this point, subject can stand up from chair.

5. To release the locked weight, press the  key again to return to the device to normal mode.

Note: Hold function can be activated before or after subject sits in chair. However, if subject finds it difficult to sit still, we recommend activating Hold after subject is seated.

C. BMI

- 1. In normal mode, press the key to enter BMI mode.
- 2. Display will show last recorded height. Left-most digit will flash.
- 3. Enter height using numeral keys (ex: 170 cm). Input will automatically move to next digit. Press key to re-input. Press key to manually move to next digit.
- 4. After inputting height, press  to confirm.
- 5. Proceed to weigh subject as usual. Indicator will display weight, height, and BMI.
- 6. **NOTE:** Hold function can be used at this time if weight is unstable
- 6. Press key to return to normal mode.

Category	BMI (kg/m ²)	Risk of obesity-related disease
Under	< 18.5	Low
Normal	18.5-24.9	Average
Over	24.9-29.9	Slightly Increased
Obese I	30.0-34.9	Increased
Obese II	35.0-39.9	High
Obese III	> 40	Very High

(World Health Organization adult BMI standards)

D. Tare

The tare function allows the user to deduct the weight of objects from the device's measurement result.

- 1. Place object that needs to be tared onto measurement platform.
- 2. Press key after stable symbol appears on indicator.
Display will indicate "0.00 kg".
- 3. Guide subject (plus tared object) to sit in chair. Conduct measurement.
- 4. To clear tare value, remove all objects from measurement platform, and press key.

E. Pre-Tare

The Pre-Tare function is used to subtract the known weight of a substance prior to weighing. The device can store 5 sets of pre-tare values.

Pre-tare values can be stored using two different methods: "Load Weight", or "Input Manually".

After pre-tare weights have been stored, they can be recalled by holding the  key for 3 seconds.

A. Load Weight

DESCRIPTION	EXAMPLE
Press  key after loading weight on the platform; the indicator will display blinking "m" symbol.	
Press numeral key 1 ~ 5 to assign this number with the current pre-tare weight.	
Press  key to store pre-tare weight; the indicator will make a beep sound.	

B. Input Manually

DESCRIPTION	EXAMPLE
Press  key. Left-most digit will begin blinking. If no further action is taken within 6 seconds, indicator will return to normal mode	
While digit is blinking: Enter pre-tare weight using 0~9 keys. Ex: to pre-tare 5.0 kg of weight, press 0-0-5-0. Ex: to pre-tare 13.5 kg of weight, press 0-1-3-5. Press  key to confirm the pre-tare weight.	
Indicator will display minus sign to the left of pre-tare weight value.	
To save this pre-tare weight value in memory: Press  key; the blinking "m" symbol will appear on the display.	

Press numeral key 1 ~ 5 to assign this number with the current pre-tare weight.



Press  key to store pre-tare weight; the indicator will make a beep sound.



NOTE: Pre-tare weight must be under max capacity, otherwise screen will show 0.00 after **[ENTER]** key is pressed, and the operator will have to re-input pre-tare settings.


C. Recall Pre-Tare Weight

DESCRIPTION	EXAMPLE
Press and hold  key for 3 seconds. Indicator will display pre-tare value M1 first. The pre-tare value will flash.	
Press numeral keys 1 ~ 5 to choose pre-tare value	
Press  key to confirm which pre-tare weight to select; the device will automatically deduct pre-tare weight.	
Press  key to return to Normal Mode	

NOTE: Pre-tare weight must be under max capacity, otherwise screen will show 0.00 after  key is pressed, and the operator will have to re-input pre-tare settings.

F. Print

If thermal printer is connected to indicator, results can be printed by pressing  key.

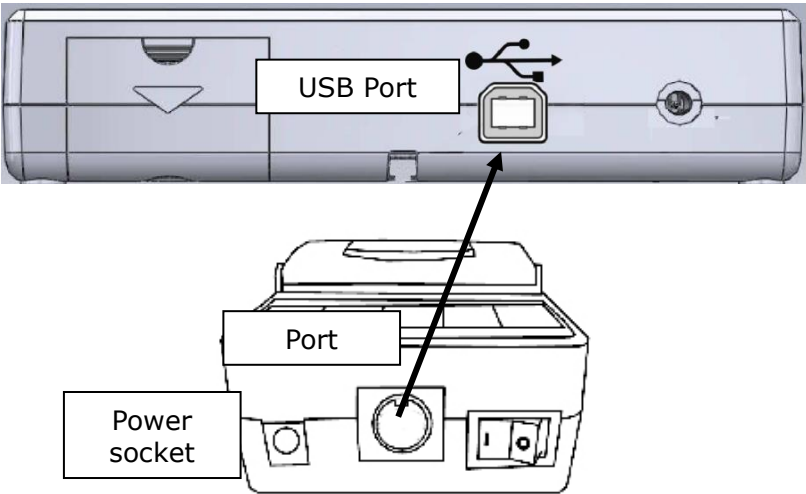

NOTE: Thermal printer needs to be powered by adapter

V. Device Setup

A. Setting Time & Date

Press and hold  key for 3 seconds to enter Time Setting mode.
Example: Inputting 2008, Dec 25, 8:00am

	Year Setting Enter year using numeral keys 0-9. Press  key once completed to proceed to month & day setting.
	Month & Day Setting. Enter month, followed by day using numeral keys 0-9. Ex: December 25th is "12.25". Input 1-2-2-5. Press  key once completed to proceed to time setting.
	Time Setting Enter time (24hr format) using numeral keys 0-9. Ex: 08:00am is input by pressing 0-8-0-0. Press  key once completed to confirm time settings and proceed to confirmation.
	Device will display new time and date settings, cycling between year, month & day, and time. YYYY→MM.DD→:HH:MM Press  key to return to normal weighing mode.

B. Device Setup

When the device is switched on, press and hold the  key for about 3 seconds, until the display shows the "SETUP", followed by "A.OFF" (first option in setting menu).

In device setup menu:



to toggle next menu option



to toggle previous menu option



to confirm selection / enter submenu



Auto Power-Off: Instruct device to shut off automatically after a certain period of time.



Press  to toggle between options (120 sec / 180 sec / 240 sec / 300 sec / off), and  to confirm selection.



Buzzer/Beep:

When function is turned on, beeping noise will be made when: indicator is turned on, keys are pressed, and weight is stable.



Press  to toggle between on/off, and  key to confirm selection.



Hold Stop: When Hold Stop is "on", Hold will deactivate after subject leaves measurement platform.



Press  to toggle between on/off, and  key to confirm selection.

Language: Set thermal printer language

Press  to toggle between English, Italian and Polish. Press  key to confirm selection.

Font size: Set thermal printer font size.

Press  to toggle between normal and double (larger). Press  key to confirm selection.

BT / Wifi (optional): If device has BT or Wifi module installed, the function can be turned OFF/BT/Wifi.

Press **[HOLD]** to toggle between OFF/BT/Wifi, and **[TARE]** to confirm selection.

Print Set (optional): If device has Wi-Fi module installed, this option will appear.

Press  to toggle between on/off, and  to confirm selection.

VI. Setup USB Connection to PC

For successful connection, PC hardware connected to device must be compatible with USB 2.0 or above. Operators should select a USB cable length that is most suitable to the operating environment.

1. Charder Smart Data Manager can be used to connect the device to a PC. The software program can be downloaded from the Charder website:

[LINK URL]<https://www.chardermedical.com/download.htm>

2. Connect USB cable to device indicator and PC. Follow installation instructions.

Program Setup

1. After installation of Charder Smart Data Manager is complete, software will automatically search for COM port. Press [**Connect**]. Once connected, [**Connect**] button will change to [**Disconnect**].

 Smart Data Manager

COM 

Connect



Gross Weight	0.0	kg	First Name	Enter
Tare Weight	0.0	kg	Last Name	Enter
Net Weight	0.0	kg	Patient ID	Enter
Height	0.0	cm	Date of Birth	31 / 12 / 1990 
BMI	0.0		Gender	Male Female
Data	Auto	Manual		

Please press "Connect".
Update Time:
Model:

Collect

Clear

Save as






Conducting Measurement

1. Input subject's first name, last name, patient ID, date of birth (DD/MM/YYYY), gender, and height (for BMI calculation) into software if needed. Press **[Clear]** to clear all input.

NOTE: information can also be input after weight measurement.

chorderSmart Data ManagerCOMConnect

Gross Weight0.0kg

Tare Weight0.0kg

Net Weight0.0kg

Height167.0cm

BMI0.0

DataAutoManual

First NameJane

Last NameDoe

Patient ID20190201

Date of Birth31 / 12 / 1965

GenderMaleFemale

Please press "Connect".

Update Time:

Model:

CollectClearSave as

2. Conduct measurement. If **[Auto]** is selected, results will be transmitted from device to software automatically and displayed on the left of screen. If **[Manual]** is selected, user must press "Collect".

chorderSmart Data ManagerCOM5Disconnect

Gross Weight72.5kg

Tare Weight0.0kg

Net Weight72.5kg

Height167.0cm

BMI26.0

DataAutoManual

First NameJane

Last NameDoe

Patient ID20190201

Date of Birth31 / 12 / 1965

GenderMaleFemale

Data updated.

Update Time: 06/03/2020 11:40:05

Model:

CollectClearSave as

Saving & Printing Results

1. Press **[Save as]** to save measurement results as .csv file on PC. Default file name is same as user ID. (ex: 20190201.csv) To track changes and multiple measurements for the same subject, we recommend not changing the default file name.

 Smart Data Manager COM 5 Disconnect

Gross Weight72.5kg

Tare Weight0.0kg

Net Weight72.5kg

Height167.0cm

BMI26.0

DataAutoManual

First NameJane

Last NameDoe

Patient ID20190201

Date of Birth31 / 12 / 1965

GenderMaleFemale

 Data updated.
Update Time: 06/03/2020 11:40:05
Model:

Collect

Clear

Save as



2. Result example:


	A	B	C	D	E	F	G	H	I	J
1	Patient ID	First Name	Last Name	Date of Bi	Gender	Gross Weig	Tare Weight	Net Weight	Height	BMI
2	20190201	Jane	Doe	31/12/1965	Male	72.4 kg	0.0 kg	72.4 kg	167.0 cm	26
3										
4										
5										

If previous results were saved in "20190201.csv", new results also need to be saved as "20190201.csv" (overwriting old file) in order to save multiple results for the same subject.

	A	B	C	D	E	F	G	H	I	J
1	Patient ID	First Name	Last Name	Date of Bi	Gender	Gross Weig	Tare Weight	Net Weight	Height	BMI
2	20190201	Jane	Doe	31/12/1965	Male	72.4 kg	0.0 kg	72.4 kg	167.0 cm	26
3	20190201	Jane	Doe	31/12/1965	Male	75.2 kg	0.0 kg	75.2 kg	167.0 cm	27
4										

Results will be saved in chronological order of measurement.

3. Press the printer icon to print out result using a printer connected to the PC.

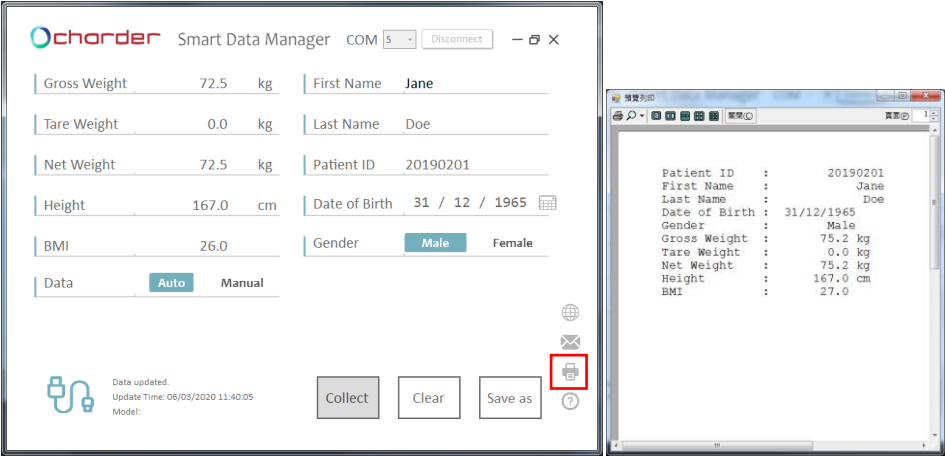

VII. Wireless Connection

If the device has the wireless or bluetooth module installed, the indicator can transmit measurement results wirelessly. Please see Charder wireless or bluetooth software instructions for details.

VIII. Troubleshooting

Before contacting your local Charder distributor for repair service, we recommend considering the following troubleshooting procedures:

Self-inspection

1. Device will not power on

- If battery power is depleted, replace with new batteries
- If batteries are not used, check if the power adapter is plugged into the device properly. Check if power adapter is plugged into mains properly

2. Indicator showing "0000" ZERO SPAN out of range

- Interference due to factors such as RF disturbance or ground vibration. Relocate device to location without interference and try again
- Unstable platform feet - adjust wheel level according to bubble level indication and try again
- External objects interfering with measurement platform. Clear platform of objects and try again
- Device may not function properly on soft surfaces such as carpets or lawns. Relocate device to location with solid, stable floor
- If the steps above cannot resolve the problem, re-calibration may be required to correct weighing accuracy

3. Connection failure for data transmission to PC or printer

- Ensure wires are connected correctly between indicator and PC or printer
- Ensure printer is supplied with power. Ensure PC software is set up properly as indicated in this manual

Distributor support required

If the following errors occur, we recommend contacting your local Charder distributor for repair or replacement services:

1. Device will not power on

- Faulty on/off key
- Broken or damaged wires causing short circuit or faulty connection

- Safety fuse burnout
- Faulty Adapter

2. Indicator damage

- Possible hardware defects include: uneven brightness in LCD screen, blurred text, smeared rainbow screen, incorrect decimal display
- Unable to save or read data
- Indicator shows "ERRL" after device is switched on
- Keys not responding
- Buzzer malfunction

Error Messages

Error Message	Reason	Action
	Low battery warning Voltage of battery is too low to operate device	Replace batteries, or plug in AC adapter
	Overload Total load exceeds device's maximum capacity	Reduce weight on measurement platform and try again
	Counting Error (too high) Signal from loadcells too high	Error normally caused by faulty loadcell or wiring. Please contact distributor
	Counting Error (too low) Signal from loadcells too low	Error normally caused by faulty loadcell or wiring. Please contact distributor
	Zero count over calibration zero range +10% while power on	Re-calibration required. Please contact distributor
	Zero count under calibration zero range -10% while power on	Re-calibration required. Please contact distributor
	Program Error Fault with device software	Error normally caused by faulty loadcell or wiring. Please contact distributor

IX. Product Specifications

A. Device Information


Model		MS5410	
Display		DP3710	
Weight Measurement	Capacity	300kg x 0.1kg,	
	Accuracy	±1.5e	
	OIML	OIML- Class III	NON-OIML
	LCD Screen	1.0-inch LCD screen (5 1/2 digits)	
Dimensions	Overall	950(W) x 1040(D) x 640(H) mm	
	Seat	Height: 560mm Width: 440 mm Depth: 360 mm	
	Armrest	Height: 700 mm	
Device Weight		22.4 kg	
Key Functions		On/Off, Zero, Print, BMI, Hold, Pre-Tare, Tare, Clear, Enter, 0~9, M1-5	
Data Transmission		USB, Wireless (optional) NOTE : Device should be connected to network by qualified distributors only.	
Power Supply		Rechargeable battery pack (optional) or 6 AA batteries / Power adapter	
Operation Environment		0°C~+40°C 15% / 85% RH 700 hPa ~1060 hPa	
Standard Accessories		User manual x 1 Power Adapter x 1 USB transfer cable x1	
Optional Accessories		Thermal Printer	

B. Power Adapter Standards



Warning

The device is only compatible with the power adapters listed below.

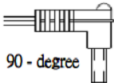AMP VOLTAGE	DRAWING NO.	CE APPROVED TYPE NO. / MODEL NO.	TY PE	Adapter plug
12V 2A	CD-AD-00041	UES24LCP-120200SPA	US	 90 - degree
	CD-AD-00041	UES24LCP-120200SPA	EU	
	CD-AD-00041	UES24LCP-120200SPA	UK	
	CD-AD-00041	UES24LCP-120200SPA	AU	

Notes

This image shows a single sheet of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.

Notes

[illegible]

X. Declaration of Conformity

Manufacturer hereby declares that this product is in conformity with the regulations and standards outlined in the following directives:

	(EU) 2017/745 Regulation on Medical Devices
	2014/31/EU Non-automatic Weighing Instruments Directive

RoHS Directive 2011/65/EU and Delegated Directive (EU) 2015/863

Radio equipment and telecommunications terminal equipment Directive 2014/53/EU
(applicable if device has wireless function)

Part 15 of the Federal Communications Statement Rules

This device may not cause harmful interference
This device must accept interference received without causing undesired operation

(Full Declaration of Conformity is available on manufacturer’s website)

Authorized EU Representative:



Obelis s.a.
Bd Général Wahis, 53
B-1030 Brussels
Belgium

Manufactured by:



Charder Electronic Co., Ltd.
No.103, Guozhong Rd., Dali Dist.,
Taichung City 41262 ,Taiwan

CD-IN-00384 REV008 12/2025