

HM110M


Infant Measuring Mat

Measuring baby length is as easy as laying them down on top of the HM110M Infant Measuring Mat. The soft, flexible, vinyl mat is printed on both edges with measurements available in centimeters and inches. The rigid footboard folds to the bottom of the infant's feet, and works with the headboard to ensure accurate measurements, with instructions printed directly on the mat.

Professionals will appreciate that the Infant Measuring Mat is portable, lightweight, and easily stored.

Features

- Graduation in inches and centimeters along both edges for easy viewing
- Non-toxic and safe for baby skin
- Can be cleaned with medical detergent
- High-quality screen print
- Instructional diagram printed directly on mat

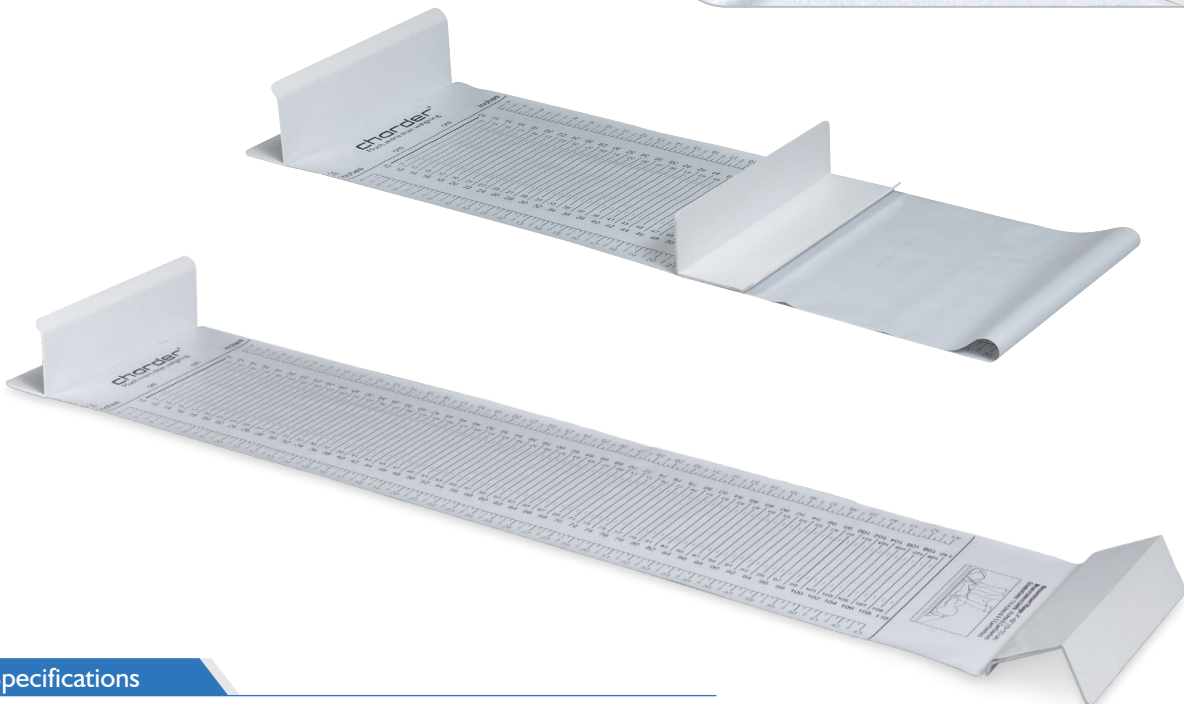

Specifications

Capacity	10-110 cm (3.9-43.3 inch)
Graduation	5mm (1/4in)
Dimensions	1400 (L) x 300 (W) mm
Device Weight	0.4 kg