

# MHS2710

## Lift Scale

A cost-effective method for weighing patients with limited mobility, this compact lift scale fits neatly between the lift system and sling, without significantly affecting the height of the lift.

Featuring a rechargeable lithium battery pack and wireless data transmission for improved workflow, the MHS2710 is a completely self-contained lift scale designed to be operated quickly and easily.



### High Capacity Rechargeable Battery

Rechargeable battery provides up to 160 hour operating time



### Battery Quick-Swap

Quick-swap system allows battery to be swapped out, minimizing downtime



### Wide Viewing Display

Featuring a new wide angle LCD indicator, measurement results can be viewed clearly from all angles




### Wireless Transmission

With wireless module installed, results can be sent directly to EMR system to simplify workflow




### Versatile Compatibility

Several joint attachments are available, making the Lift Scale compatible with a wide variety of systems

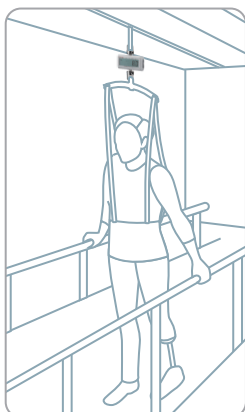

Versatile usage

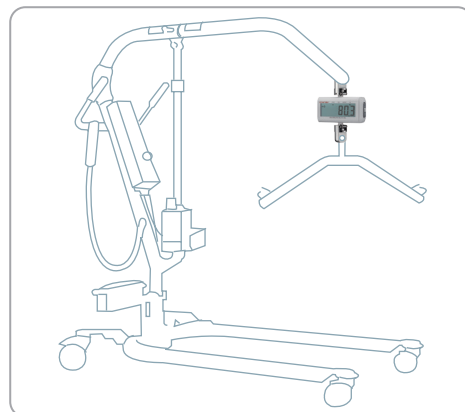

Compatible with different life systems



Multiple joint connectors




## Part of an Integrated System

After completing weight measurement, results can be sent wirelessly to a tablet, for direct record keeping or further export into EMR systems if needed.

## Features



- BMI
- HOLD
- TARE
- Battery Quick-swap
- Real-time result sync
- Up to 160 hours usage
- Power level indicator
- Wireless transmission
- Data can be managed using Charder Connect app

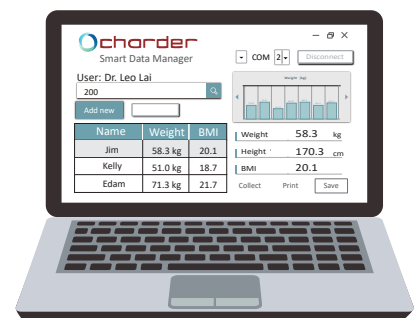

Wireless

- Weight
- Height
- BMI



Capacity / Graduation	300kg x100g / 400kg x 200g
Key Functions	On/Off/Tare/Zero, Hold/BMI
Screen Dimensions	95 x 34 mm
Overall Dimensions	122(W) x 60(D) x 180(H) mm
Device Weight	1.04 kg
Power Source	AAA x 4 Rechargeable Battery (optional)
Accessories	Micro USB charging cable (optional) Rechargeable Battery (optional)



Charder Electronic Co., Ltd.  
 No.103, Guozhong Rd., Dali Dist., Taichung City 41262 Taiwan  
 TEL: +886 4 2406 3766 FAX: +886 4 2406 5612  
 Email: info\_cec@charder.com.tw www.chardermedical.com