

MS3830

Wheelchair Scale

The high-capacity MS3830 wheelchairscale is portable and easily transported on rolling castors. Optional handrails can be attached for added stability when used by patients who have difficulty standing. By deducting the weight of the wheelchair using the scale's tare function, a patient's exact weight can be accurately measured.



Image contains optional handrail



Castor wheels for portability

Functions


- BMI
- Tare
- Portable Handle
- Print
- Pre-Tare
- Wheel Castors
- Hold
- Backlight
- Data Transfer

Specifications

Capacity	300 kg
Graduation	0.1 kg
Overall Dimensions	1150(W) x 800(D) x 66(H) mm (not including handrails)
Platform Dimensions	900(D) x 740(D) mm
Device Weight	28.6 kg
Power Supply	AA batteries*6 / Power Adapter / Rechargeable battery pack (optional)
Optional Accessories	Thermal Printer, Display Stand, Handrails



Push wheelchair onto platform
via two-way ramps

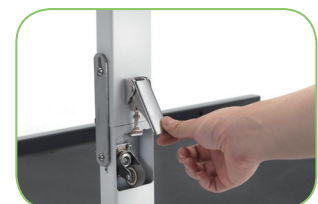

Quick lock and release for
handrails