



# MS4640


## Portable Medical Scale with remote display

The MS4640's remote display can be handheld or mounted to an optional stand or the wall, and records weight, height (manually entered), and BMI. USB connectivity allows data transmission to a PC or an optional thermal printer. Wireless modules are also available.

This scale features a large, durable, die-cast metal platform with an anti-slip rubber mat for safe and accurate weighing. A bubble level and height-adjustable feet allow users to make the platform of the scale even for patients to stand on safely.



Portable Handle

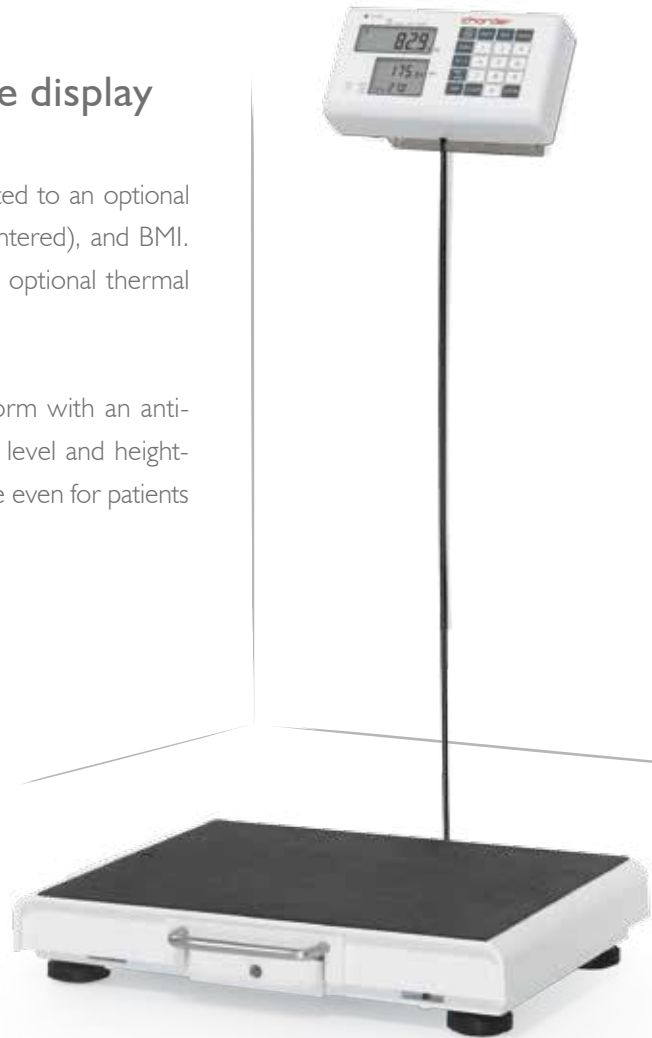

### Functions

- BMI
- Print
- Hold
- Backlight
- Tare
- Pre-Tare
- Portable Handle
- Remote Display
- Bubble Level
- Adjustable Feet
- Anti-Slip Mat

### Specifications

Capacity	300 kg
Graduation	0.1 kg
Overall Dimensions	450(D) x 340(W) x 90(H) mm
Device Weight	8.4 kg
Accessories	AR-249I Carrying Bag TP2100 Thermal Printer
Data Transfer	USB, Wireless
Software	Smart Data Manager, Charder Connect



AR-249I Carrying Bag